



ACTGene, Inc.
PO.Box 7011
Piscataway, NJ 08855-7011
Tel: 973-978-7888
Email: sales@actgene.com
Web Site: www.actgene.com

ACTGene™ Protein Marker 10-180 kD, 10 Protein 3 Color

Product contents:

P/N: ACT-IDMW26

1x500µl

Sufficient for approximately 150 loads

Introduction:

ACTGene™ Protein Marker is ten protein three color pre-stained marker covering wide range of 10-180kD protein sizes. Conveniently delivered in ready-to-use format and ideal for SDS-PAGE, Western Blotting, and related protein applications.

Protocol:


Swirl tube well.

Add ~4µl per load onto protein gel prior to electrophoresis.

Shipping and Storage Conditions:

Shipped in RT

Store in in 4°C for short term and -20 °C for long term

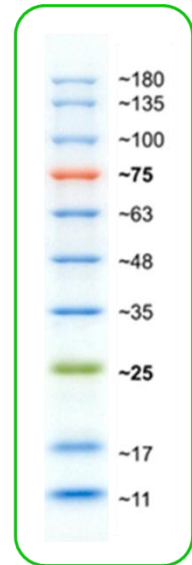

Tris-Glycine 4-20%
(10-180kD) 10 Protein 3 Color