

## ACTGene™ DNA Marker 100 bp Ladder, 11 Fragments

### Product contents:

P/N: ACT-IDMWD100

1x500µl

Sufficient for approximately 100 loads

---

### Introduction:

ACTGene™ DNA Marker is eleven DNA fragment ladder covering wide range of 10-1,500 bp sizes. Conveniently delivered in ready-to-use format and ideal for agarose gel electrophoresis, Northern, Southern, and related DNA analysis applications.

---

### Protocol:

Swirl tube well.


Add ~5µl per load onto agarose gel prior to electrophoresis.

---

### Shipping and Storage Conditions:

Shipped in RT

Store in in 4°C for short term and -20 °C for long term

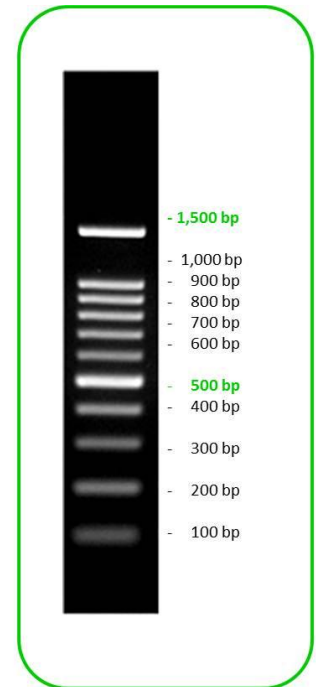

1.5% Agarose Gel 1xTAE  
(100-1,500 bp) 11 Fragments